

i i i i i i i i i i i i i i i i i i i	Face	Macr
ı un	Iacc	i'iasr

Might be right for you if you	Are a mouth breather
	Experience chronic nasal obstruction or congestion
Might not be right for you if you	Feel claustrophobic





# Fitting guide





Step 1.

Hold the front of the mask with one hand and the headgear with the other. Ensure that one headgear clip is unhooked from the frame.



Step 3.

Hook the unattached headgear clip onto the frame.



Step 2.

Place the seal onto your face and guide the headgear over your head.



Step 4.

Gently tighten the headgear straps of the headgear starting with the bottom straps, followed by loosely fitting the blue forehead straps.

# Reorder codes

### Masks


Code	Product
VIT1SU	Vitera Mask Full Face Small
VIT1MU	Vitera Mask Full Face Medium
VIT1LU	Vitera Mask Full Face Large

Code	Product
400VIT111	Vitera Seal Small
400VIT112	Vitera Seal Medium
400VIT113	Vitera Seal Large
400VIT121	Vitera Headgear Small
400VIT122	Vitera Headgear Medium/Large
400FPH241	F&P Swivel
400VIT151	Vitera Clips and Forehead Clip
400VIT171	Vitera Elbow







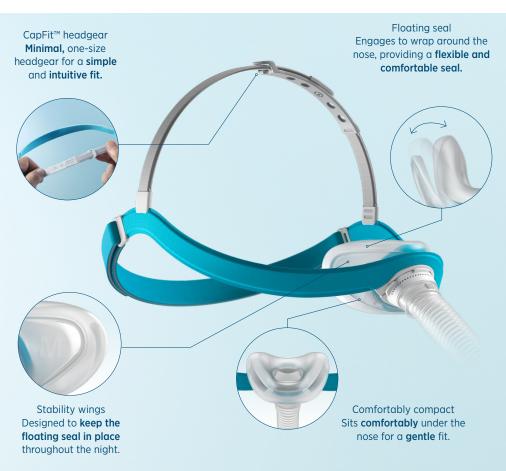

Compact Nasal Mask
Don't want something in your nose
Want a clear line of sight
Need a mask that is easy to fit
Want a smaller mask
Are a mouth breather











For further instructions, please refer to Fitting Your Mask section of the Evora Use and Care Guide.

### Reorder codes


Masks	
Code	Product
EVO1SA	Evora <sup>TM</sup> Nasal Mask Small A Model
EVO1MA	Evora Nasal Mask Medium A Model
EVO1LA	Evora Nasal Mask Large A Model
EVO1WA	Evora Nasal Mask Wide A Model
EVO1SMLU	Evora Nasal Mask Fit Pack SML

Code	Product
400EVO115	Evora Nasal Seal Small Spare
400EVO116	Evora Nasal Seal Medium Spare
400EVO117	Evora Nasal Seal Large Spare
400EVO118	Evora Nasal Seal Wide Spare
400FPH341	F&P Swivel White
400FPH181	F&P Tube Strap
400EVO121	Evora Nasal Headgear Spare
400EVO122	Evora Nasal Backstrap Spare







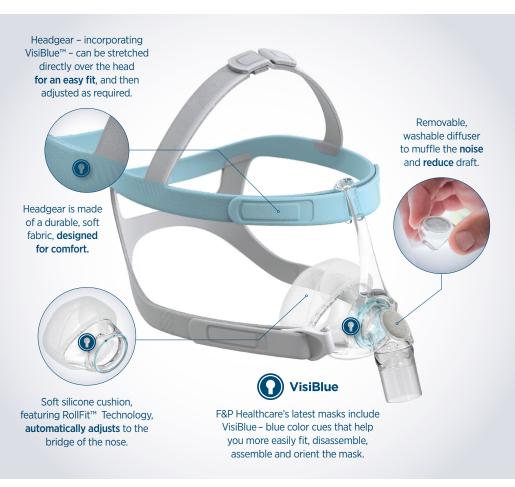

Minks by Minks Committee of the Committe	Don't want something in your nose
Might be right for you if you	Need a mask that is easy to seal and fit
Might not be right for you if you	Are a mouth breather





# Fitting guide





Step 1.

Hold open the bottom straps of the headgear with the mask frame facing downwards.



Step 3.

Place your hands on the bottom straps on either side of the mask frame and pull out and downwards. Ensure the silicone seal rests over your nose and the blue top strap is sitting on your forehead.



Step 2.

Guide the headgear over your head so that the silicone seal is resting on your forehead.



Step 4.

Gently tighten the headgear straps of the mask headgear starting with the blue forehead straps followed by the bottom straps.

### Reorder codes

#### Masks


Code	Product
ESN2SA	ESON2 Small
ESN2MA	ESON2 Medium
ESN2LA	ESON2 Large

Code	Product
400ESN214	ESON2 Seal Small
400ESN215	ESON2 Seal Medium
400ESN216	ESON2 Seal Large
400ESN221	ESON2 Headgear Small
400ESN222	ESON2 Headgear Medium-Large
400FPH241	F&P Swivel 2
400ESN251	ESON2 Clips
400ESN261	ESON2 Diffuser
400ESN271	ESON2 Elbow







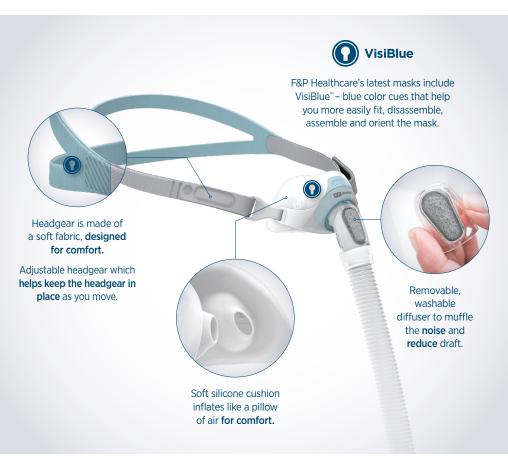

	Nasal Pillows Mask
	Feel claustrophobic
Might be right for you if you	Want a clear line of sight
	Want a smaller mask
Might got be wight for you if you	Have very sensitive nostrils
Might not be right for you if you	Are a mouth breather





# Fitting guide





Step 1.

Hold the silicone seal away from your nose and use your fingers to carefully spread open the soft side supports of the silicone seal and then guide the two openings into your nostrils.



Step 3.

Gently tighten or loosen the headgear straps on the sides of the headgear until you achieve a comfortable seal.



Step 2.

Hold the mask frame in position with one hand and pull the headgear over your head with the other hand, positioning the straps above the ears as per the fitting image above.

# Reorder codes

#### Masks

Code	Product
BRE1SA	Brevida Extra Small-Small (one XS-S seal)
BRE1MA	Brevida Medium-Large (one M-L seal)
BRE1SMA	Brevida XS-S and M-L (two seals)

Code	Product
400BRE113	Brevida Seal Extra Small-Small
400BRE114	Brevida Seal Medium-Large
400BRE120	Brevida Headgear
400BRE151	Brevida Headgear Clips
400BRE161	Brevida Diffuser
400BRE171	Brevida Elbow
400FPH241	F&P Swivel 2

